

Powered
by



PERFECT
DENT

KING RAPTOR ENDOMOTOR



4

YEARS
WARRANTY


IP MAKING
THE FUTURE

RED IS TAKING
OVER EGYPT

BRUSHLESS MOTOR

ADVANCED MAGNETIC TECHNOLOGY

- Reduce mechanical friction
- Quiet less noise compared to brushed motors
- Longest life span ever
- Highest efficiency & performance
- Consumes less power and provides consistent torque
- **Lower heat production:** comfortable use despite extended working hours
- High precision ensures accurate control of speed and torque

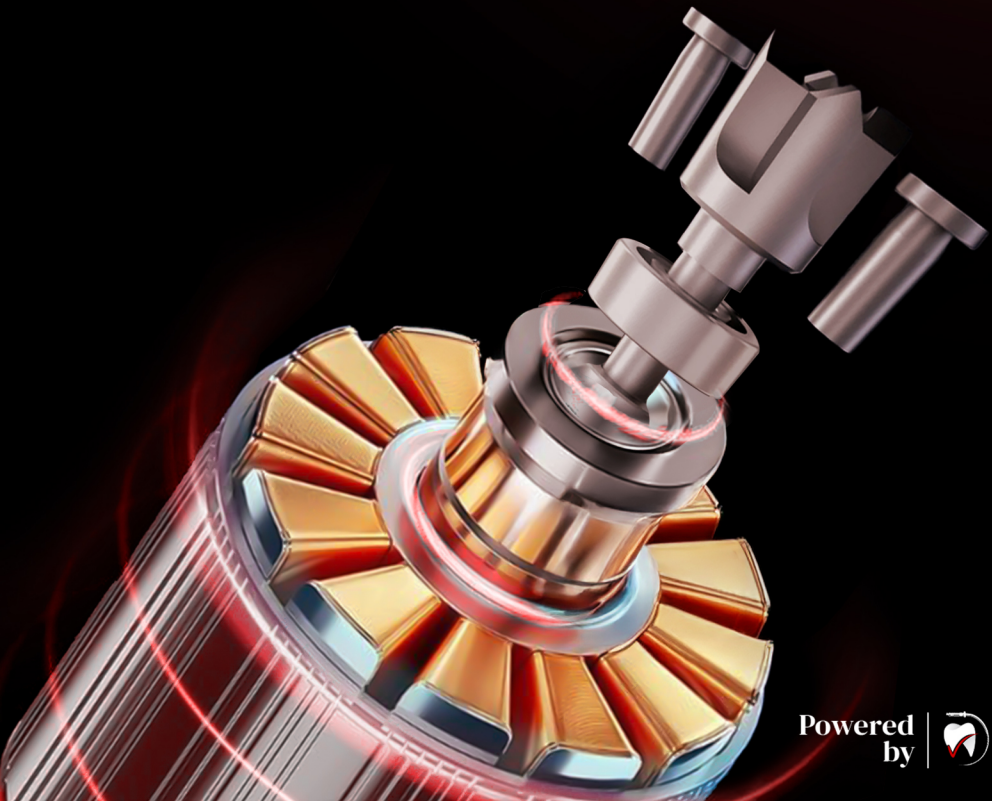

LUXURY DESIGN, MASTER PERFORMANCE



MINI HEAD CONTRA, SLIM DESIGN



- **Reinforced Stainless Steel Contra:**
 1. Maximum rust and corrosion resistance
 2. Maximum impact resistance
- **6:1** Gear reduction provides consistent torque and speed whatever the pressure
- **360°** rotatable contra angle

PREMIUM ERGONOMIC DESIGN

- **Weight:** well balanced only 160g
- Design for maximum ergonomic performance
- Richness shiny body reinforced with metal ornament
- **OLED Unmatched Screen:**
 1. Premium futuristic viewing edges provides visibility from multiple angles for more immersive experience.
 2. Lighter screen with high resolution regardless of room/unit lighting
 3. Energy efficiency consumes less power
 4. Slimmer elegant display





Powerful Battery: 1900 mAh



Speed: 100 - 2000 rpm



10 programs



Torque: 0.5 - 6 N.cm

L/R

Suitable for right & left handed operators



Motions:

1. CW (clockwise)
2. CCW (counter clockwise)

3. ATC (clockwise):

- * Fully opened adjustment of angles (10° - 360°)
- * Speed starts from 100 rpm
- * Torque starts from 0.6 N.cm

4. REC (reciprocation)

- * Fully opened adjustment of angles (10° - 360°)
- * Speed starts from 100-600 rpm
- * Torque starts from 2.0 - 6.0 N.cm

Safety Features:



Auto Reverse



Adaptive torque control



Safety Features: Prevents fracture - Sensitive faster response



New buttons withstand countless presses



MICRO USB type C, rechargeable, fast charging


CONTACT US



perfectdenteg.com



perfectdenteg

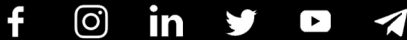

perfectdenteg



info@perfectdenteg.com



+201552974571

+201018436971



74 Ammar ibn Yasser, above Metro Market,
Masr Al Jadidah, Cairo, Egypt